

iProphy COMPLETE SETS



MODEL : **iProphy Complete set**
ORDER CODE : **Y1500015**

- Contents :
- iProphy Motor handpiece • iFX57 Head
 - iProphy Stand • AC ADAPTOR • F-Type Spray Nozzle

iProphy without Head

MODEL : **iProphy W/O**
ORDER CODE : **Y1500069**

- Contents :
- iProphy Motor handpiece
 - iProphy Stand • AC ADAPTOR • F-Type Spray Nozzle

iProphy HEADS



MODEL : **iFX57**
ORDER CODE : **C1102**

- For Prophy • For screw-in type
- 4:1 Reduction



MODEL : **PR-F**
ORDER CODE : **H209**

- For disposable prophy angle
- 4:1 Reduction



SPARE PARTS

MODEL	PRODUCTS	ORDER CODE
iProphy Motor handpiece		EA11170020
iProphy Stand		UA10910020
AC ADAPTOR		Y1500014
BATTERY PACK		E1080052
F-Type Spray Nozzle		Z256090



Washable in the thermodisinfectant



Autoclavable



iProphy

CORDLESS POLISHING HANDPIECE

Put The Control Back in Your Hands.

With no foot control, the cordless iProphy from NSK gives you maximum freedom and stability for smoother polishing and greater control. Take your clinical practice to the next level by improving your performance day after day.


The cordless handpiece guarantees total freedom of movement.

Working cordless means no restrictions on maintaining a good posture and access angles resulting in less stress for operator and patient. Total freedom to follow the user to any location and treatment position.



An ultra mini head guarantees the best possible field of vision.

The thin, long ultra mini head ensures a wide field of vision within the oral cavity, easy access and excellent reach. Facilitates even smoother and more precise treatment.

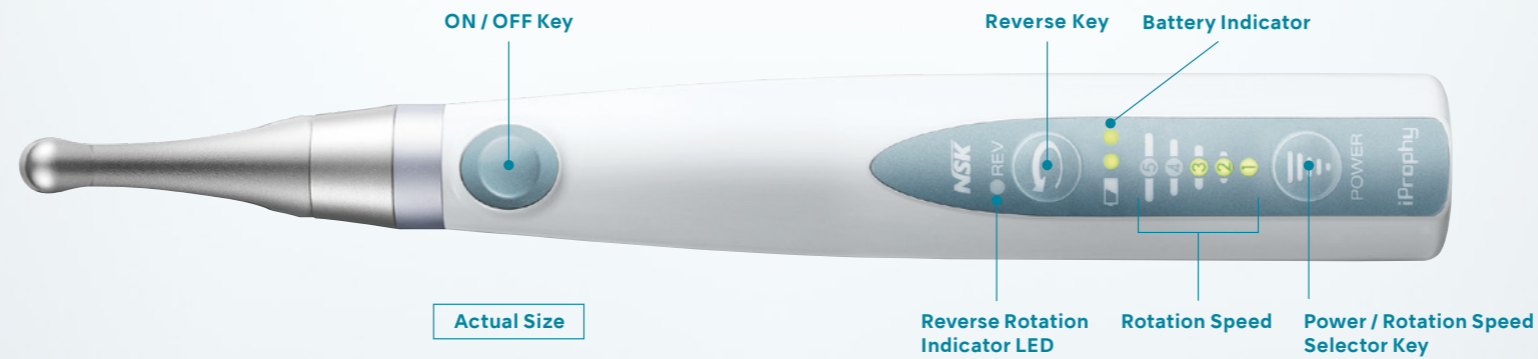

Easy power control.

Power control at your fingertips. No foot-control required. iProphy easily maintains low speed rotation, which can be difficult with conventional contra-angles.



Lightweight and easy to hold.

Ergonomically-designed and developed to include the latest ergonomic features. The lightweight, easy to hold and well-balanced handpiece feels light in the hand even when used for a long time.



Simple head switching.

Simple head switching using the quick joint system. In addition, operability is vastly improved because the head can be positioned in six different angles to best suite different procedures and easily access all areas of the oral cavity.



It has a simple and clearly designed user interface.

The intuitive user interface makes procedures stress free and contributes to a relaxed atmosphere. The ON/OFF key is large and easy to operate.



iProphy offers a wide-range of speed settings.

It is possible to set a wide-range of speeds from 500-2,500 min⁻¹. iProphy can be used for a wide-range of procedures.